


## T2 Ratio Plug-in

Antoine Rosset, 2004

Pre and Post contrast T2 ratio


1) Open two series:

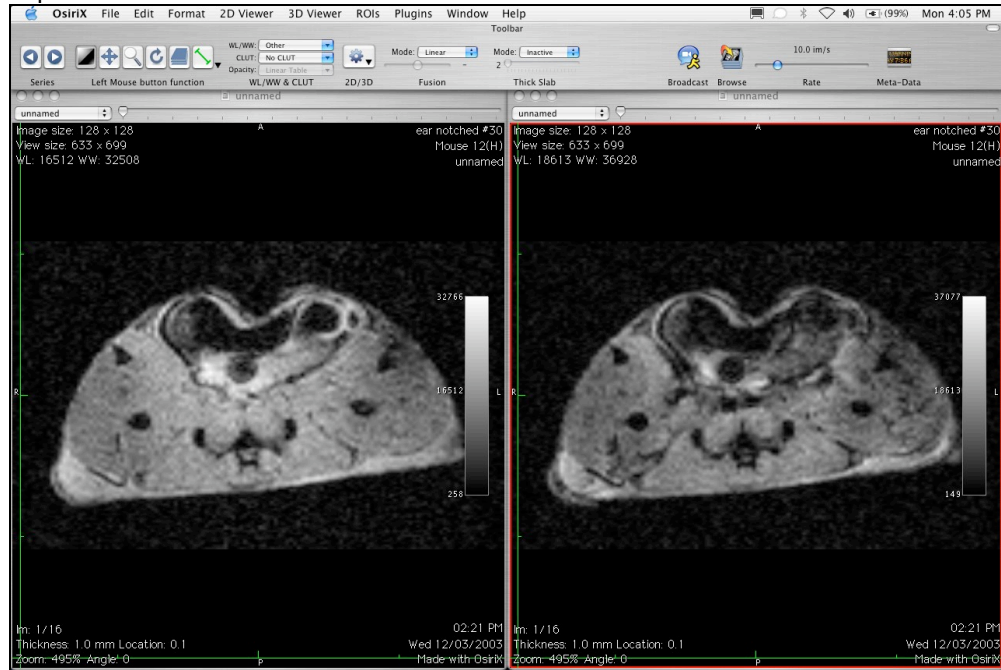


Left is PRE-contrast GRE

Right is POST-contrast GRE

2) Drag & Drop POST-contrast on PRE-contrast series. Select the **T2 Ratio** plugin.

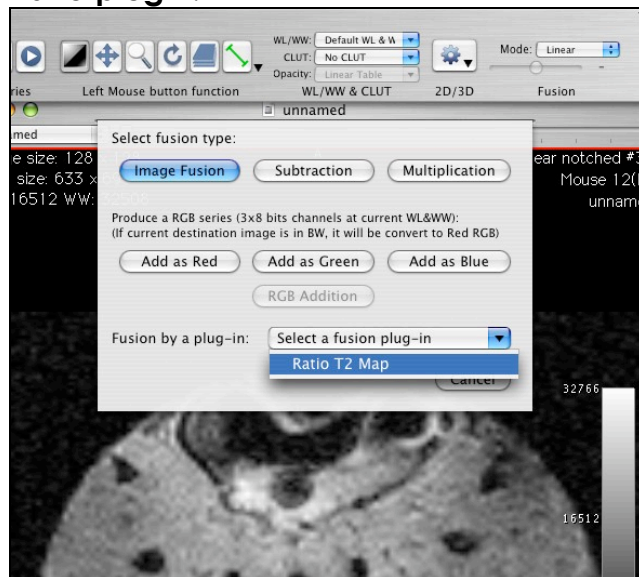

- 3) Adjust parameters in T2 Ratio window. You can limit the computation to only selected ROIs of the PRE-contrast series.

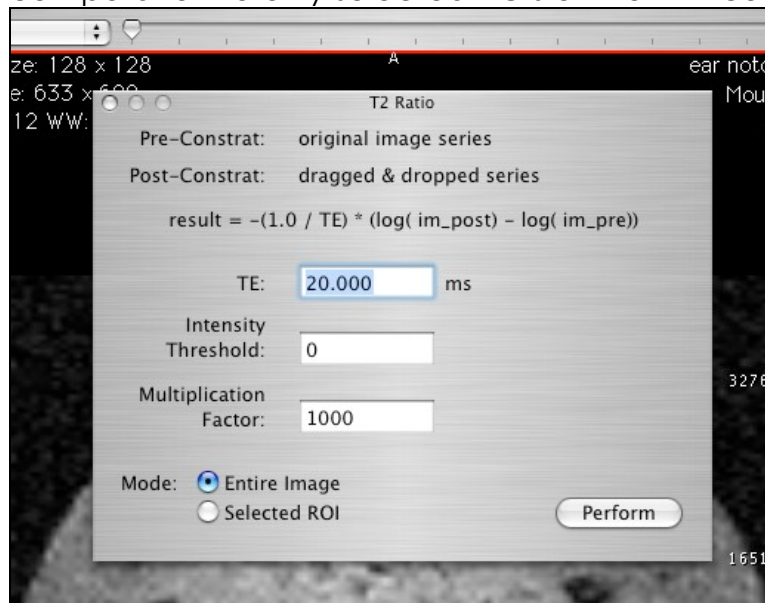

Click **Perform** to produce the T2 Ratio series.

- 4) Drag & Drop this new series on your PRE or POST contrast series to produce a series with image fusion. Adjust the CLUT, window levels and image fusion percentage.

